suavinēx

Manual de instrucciones

Instruction manual

Manuel d'instructions



Manuale d'istruzioni

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**APP
CONTROL**

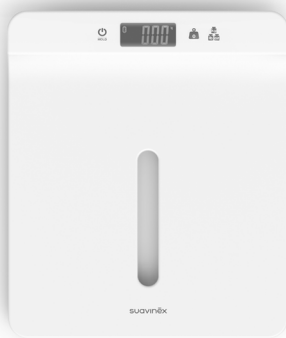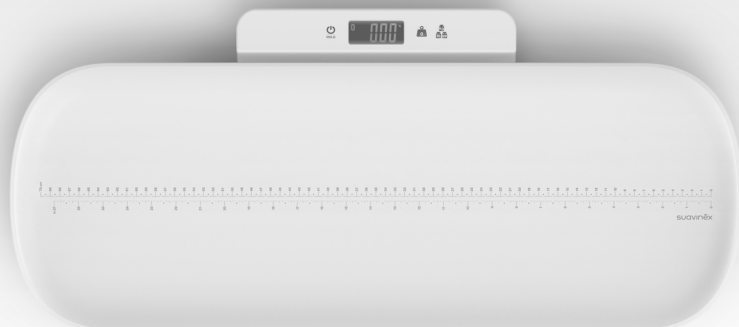ES	ADVERTENCIAS	4 - 6
EN	WARNINGS	6 - 8
FR	AVERTISSEMENTS	8 - 10
IT	AVVERTENZE	10 - 12
DE	HINWEISE	13 - 15
PT	ADVERTÊNCIAS	15 - 17

ES	MANUAL DE INSTRUCCIONES	20 - 25
EN	USER MANUAL	26 - 31
FR	MANUEL D'INSTRUCTIONS	32 - 37
IT	MANUALE D'ISTRUZIONI	38 - 43
DE	GEBRAUCHSANWEISUNGEN	44 - 49
PT	MANUAL DE INSTRUÇÕES	50 - 55

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EN: WARNING

Please read these instructions carefully and keep them for future reference.

1. This device is designed as a “weighing scale” for domestic use; any other use is considered unsuitable and therefore hazardous. The manufacturer takes no responsibility for any damage caused by improper, wrong or irrational use of the device.

2. This product is not a toy. Keep out of reach of children. Do not let children play with the device. It must NOT be used by children, people with impaired physical, sensory or mental capabilities or people with learning difficulties, unless they are supervised or have been properly informed about how to handle the product correctly and they are aware of the potential hazards. It must be used, cleaned and maintained by an adult.

3. Never leave the baby unattended during measuring.

4. The scale is powered by two AAA batteries.

5. The batteries must be inserted by an adult.

6. Do not short-circuit the battery terminals.

7. Place the product on a stable, flat surface.

8. To clean the scale, use a soft cleaning cloth, moistened with a little water if necessary.

9. For cleaning and maintenance of the mat, follow the washing instructions on the label.

10. Do not use aggressive or abrasive detergents or solvents. Do not use chemical solutions.

11. After removing the scale from its packaging, make sure it is intact and there is no visible damage. Keep the packaging out of the reach of children. If in doubt, do not use it and contact the manufacturer or distributor.

12. Certain basic rules must be observed when using any electrical device:

- Do not expose the device to atmospheric conditions such as rain, sun, etc. This scale is for indoor use only.

- Keep the device away from heat sources.

- Do not touch the device with wet or damp hands.

- Do not place or use the product in damp places.

- Do not knock the product and avoid dropping or damaging the product when moving it, as it contains high-precision electronic parts.

- Do not use a chemical solution to clean the product.

- Do not touch the device if it has fallen in water.

- Do not use it again until you have contacted the manufacturer, distributor or authorised repair centre.

13. This scale can measure weights between 50 g and 100 kg. Do not use it if the weighing capacity is exceeded.

14. The scale is not suitable for use in environments where flammable gases and/or anaesthetics are present, as it could cause explosions.

15. Do not throw the scale on a fire, as it could explode.

16. Avoid the scale coming into contact with any liquid.

17. If the scale is not going to be used for a long time, it is advisable to remove the batteries.

18. Use only original accessories, parts or components.
19. If the device is faulty or does not work properly, switch it off immediately. Do not try to repair it and contact the manufacturer, distributor or authorised repair centre. Any repairs must be carried out by authorised technical personnel.
20. Product characteristics may be changed without prior notice.

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EN: INSTRUCTION MANUAL

· CONTENTS

1. INTRODUCTION	26
2. PARTS OF THE SCALE	26-27
3. INSTRUCTIONS FOR USE	27-30
3.1. INSERTING BATTERIES	27
3.2. ASSEMBLING THE TRAY ON THE SCALE	28-29
3.3. USING THE SCALE (WITH OR WITHOUT THE TRAY) WITHOUT THE APP	29
3.4. DISMANTLING THE TRAY	29
3.5. USING THE SCALE (WITH OR WITHOUT THE TRAY) WITH THE APP	29
3.6. LOCK FUNCTION	29-30
3.7. TARE FUNCTION	30
4. CLEANING AND MAINTENANCE	30
5. TROUBLESHOOTING	30-31
6. TECHNICAL SPECIFICATIONS	31

1. INTRODUCTION

The Suavinex Smart evolving scale allows you to accurately monitor the baby's growth from the first days (both weight and length*) to the next stages of growth. With its ProWeighing technology, it allows you to accurately monitor weight. Designed for use from birth up to adulthood (max. 100 kg).

*Up to 70 cm/27 in.



You can also connect to the Suavinex Link® App via Bluetooth to save the weight, height/length and head circumference measurements of several users.

Discover all the application's possibilities at www.suavinex.com


2. PARTS OF THE SCALE


A. SCALE

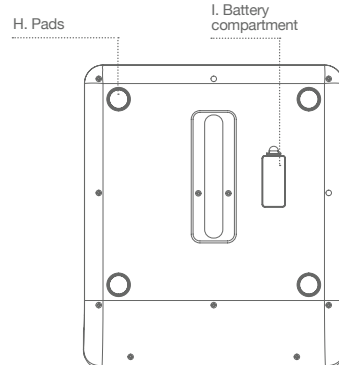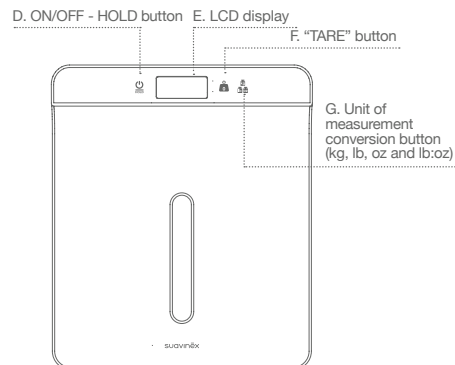B. TRAY WITH LENGTH MEASURING SCALE

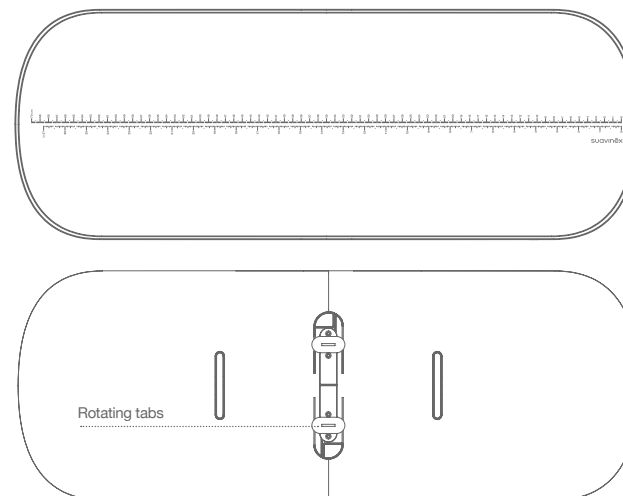C. MAT

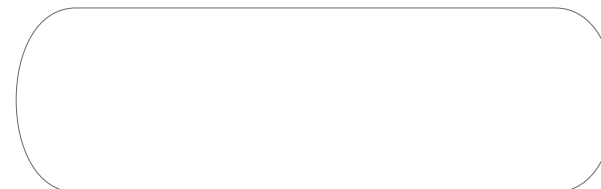
3. INSTRUCTIONS FOR USE

3.1. INSERTING BATTERIES

This scale is powered by two AAA batteries. For correct insertion, follow the steps below:
 1. Open the cover on the base of the scale.
 2. Insert the batteries, ensuring the negative and positive terminals are in the corresponding position.

-IMPORTANT:
 -When the LCD display (E) loses intensity, it means the batteries are running out. It is time to replace them.
 -If the scale is not going to be used for a long time, it is advisable to remove the batteries.
 -Place the spent batteries in a specific bin so that they can be properly recycled.

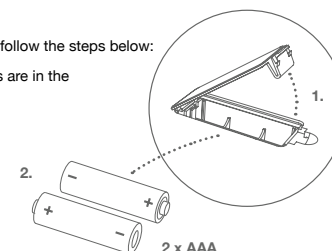3.2. ASSEMBLING THE TRAY ON THE SCALE

If the measurement is to be taken with babies, follow the steps below to attach the tray (B) to the scale (A):

1. Put the two parts of the tray together by fitting the tabs into the grooves. See FIG. 1. It is important that the rotating tabs on the back are parallel to the parts. See FIG. 2.
2. Once the two parts are assembled, carefully insert the protruding part into the slot in the base. See FIG. 3.
3. Turn the tray and scale (B+A) upside down. The rotating tabs remain in the same direction. See FIG. 3A.
4. Rotate the tabs so that they are positioned horizontally and the tray is securely attached to the scale. See FIG. 4.
5. Make sure the tray is securely attached to the scale.

IMPORTANT:

- To prevent accidents, make sure the base is securely attached.
- The tray can withstand up to 20 kg.

FIG. 1

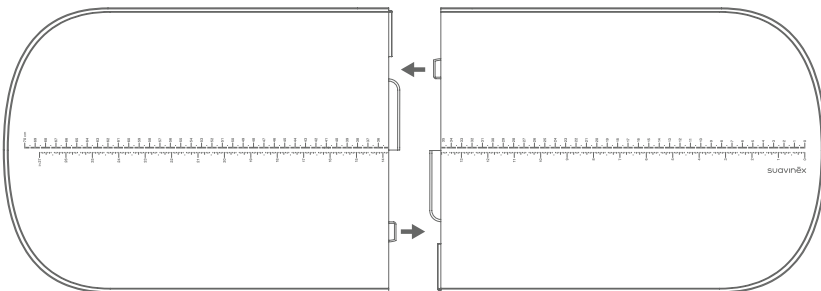
FIG. 2

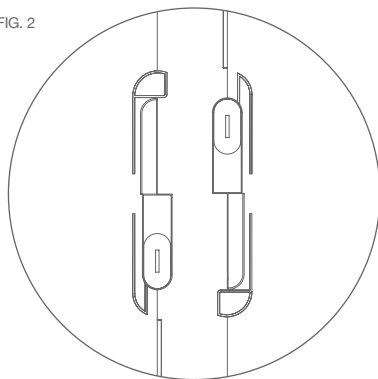
FIG. 3

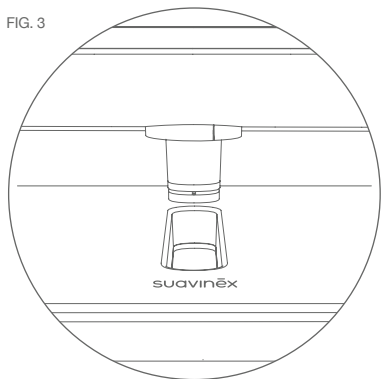
FIG. 3A


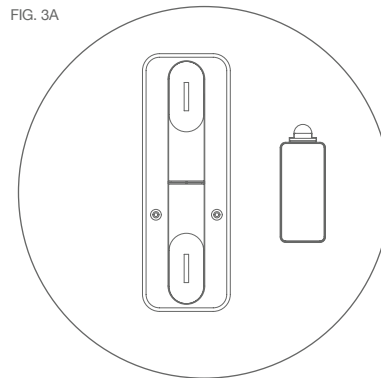
FIG. 4

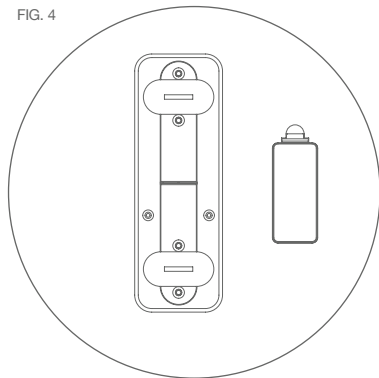3.3. USING THE SCALE (WITH OR WITHOUT THE TRAY) WITHOUT THE APP

Having attached the tray (B) or if measuring directly on the scale (A):

1. Make sure the scale is on a stable, flat surface.
 2. Press the "ON/OFF" button to switch on the scale.
 3. Once the display shows 0.00, the measurement can be taken.
 4. During measuring the weight fluctuates and once the weight has stabilised the measurement will flash and remain on the screen.
- To take another measurement, remove the weight and press the "TARE" button.

IMPORTANT:

- If the measurement is not established due to excessive movement, it is possible to establish it manually by pressing the HOLD button (1), as explained in section 3.4.
- If using the tray, it is not possible to take the measurement with the baby standing up.
- The tray can withstand up to 20 kg.
- For correct measurement, it is important for the person being weighed not to touch anything other than the scale so as not to interfere with accurate weighing.
- Up to 30 kg, the scale measures to two decimal places and the resolution is measured in 10 g increments.
- Above 30 kg, the scale measures to one decimal place and the resolution is measured in 100 g increments.
- Make sure the measurement is in the desired units of measurement (kg, oz, lb, lbxoz).
- The buttons are highly sensitive, so functions/units of measurement may be accidentally activated.
- This scale has an energy-saving mode where, if no activity is registered for 60 seconds, the scale will switch off. To switch the scale off immediately, press and hold the "ON/OFF" button until the display switches off.
- If the LCD display shows "ERR", take the baby/person off the scale and start a new measurement.

3.4. DISMANTLING THE TRAY

Follow the steps below:

- Rotate the tabs so that they are in the position shown in FIG. 3A.
- Separate the tray (B) from the scale (A) FIG. 3.

3.5. USING THE SCALE (WITH OR WITHOUT THE TRAY) WITH THE APP

In order to use the Smart scale with the Suavinex Link® App, follow the steps below:

1. Switch the scale on by pressing the ON/OFF button (A). At that point, the bluetooth icon (4) will start to flash blue. This means that the device is ready to connect to the Suavinex Link® App.
2. Enable the Bluetooth connection on your smartphone.
3. Open the Suavinex Link® App and go to the "devices" tab. All the Suavinex Smart products will appear under this tab. Select the scale.
4. Once the device is connected to the smartphone, the Bluetooth icon will remain on.

3.6. LOCK FUNCTION

The weight can be established in two ways:

- **Automatic stabilisation:** when measuring a baby or child, it is possible they may move, but thanks to the automatic locking system, when they are still enough, the weight will remain on the display. You can verify that it has stabilised when the weight flashes and remains on the display.
- **Manual stabilisation/locking:** if the baby or child moves too much, there may be problems with the automatic locking system. You can choose to do this manually by pressing the HOLD button for approximately 1 second. The

signal flashes and the weight is shown on the display. The word HOLD will be displayed on the screen. This function is available above 1 kg.

3.7. TARE FUNCTION

This function makes it possible to deduct the weight of the mat or any other textile placed on the scale:

1. Press the ON/OFF button until the scale switches on.
2. Once the scale is at "0", place the added weight on the tray (B) and the weight will be displayed on the screen.
3. Press the TARE button to restore the weight to "0".
4. The measurement can now be taken.

IMPORTANT:

-When removing the mat from the tray, the TARE button must be pressed again, otherwise a negative weight will be shown on the display and the weight obtained will not be accurate.

4. CLEANING AND MAINTENANCE

- **SCALE AND TRAY:** use a soft cleaning cloth, moistened with a little water if necessary. Do not use aggressive or abrasive detergents or solvents.
- **MAT:** follow the instructions on the label for proper maintenance.

5. TROUBLESHOOTING

PROBLEM	SOLUTION
THE SCALE DOES NOT SWITCH ON	- Check that the batteries work.
THE SCALE SWITCHES ON, DISPLAYS <i>LOW</i> AND SWITCHES OFF	- The device has run out of battery power. Replace the batteries as described in section 3.1 INSERTING BATTERIES.
THE MEASUREMENT ON THE SCALE DOES NOT STABILISE	- Make sure the scale is on a flat, stable surface. - Check all the pads are present on the base. - Make sure the baby does not move too much. - Try using the MANUAL HOLD function.
THE SCALE DOES NOT CONNECT TO THE APP	- The device may have run out of battery power. Replace the batteries as described in section 3.1 INSERTING BATTERIES. - Make sure the smartphone's Bluetooth connection is enabled. - Press the button ON/OFF button (A) for a few seconds until the Bluetooth icon begins to flash. - Follow the steps described in section 3.5. USING THE SCALE (WITH OR WITHOUT THE TRAY) WITH THE APP.
THE SCALE DISCONNECTS FROM THE SMARTPHONE	- The device may have run out of battery power. Replace the batteries as described in section 3.1 INSERTING BATTERIES. - Make sure the smartphone's Bluetooth connection is enabled. - Make sure the scale is not switched off. If it is, press the ON/OFF button (A) for a few seconds until the Bluetooth light begins to flash. - Make sure the distance between the device and your smartphone is less than 5 metres.

A NEGATIVE WEIGHT IS DISPLAYED ON THE SCREEN

-Press the TARE button as explained in section 3.5. and make sure 0.00 is displayed on the screen.

"ERR" IS DISPLAYED ON THE SCREEN

- Take the baby/person off the scale.
- Start a new measurement.

AN ANOMALOUS WEIGHT IS DISPLAYED ON THE SCREEN

- Check that the measuring units are correct (kg, oz, lb, lb:oz).

Our Customer Service department is available to answer any other questions you might have either by telephone +34 965107010 (for Spain) or by email suavinex@suavinex.com.

6. TECHNICAL SPECIFICATIONS

Standard for execution of technological parameters
product name: electronic baby scale
Model: BS01/1
Power supply: 2 Batteries *AAA
Measuring range: 50 g to 100 kg.
Resolution:
-From 50 g to 30 kg: 0.01 kg
-From 30 kg to 100 kg: 0.1 kg

Weight range	Allowed weight deviation
0-1 kg	±20 g
1-3 kg	±30 g
3-5 kg	±50 g
5-10 kg	±60 g
10-20 kg	±70 g
20-30 kg	±80 g
30-100 kg	±300 g

CE conformity:

This SCALE Model: BS01/1 complies with all the essential requirements and other provisions specified in EU Directives 2014/30/EU (electromagnetic compatibility), 2014/35/EU (electrical safety), 2011/65/EU (restriction of hazardous substances - RoHS), 2012/19/EU (waste electrical and electronic equipment - WEEE) and 2014/53/EU (radio equipment directive - RED). To request a copy of the CE declaration of conformity, please contact the manufacturer: EXCLUSIVAS RIMAR, S.L. - C/ del Rubio parcelas R83 - R94 Polígono Las Atalayas 03114 - ALICANTE - ESPAÑA/SPAIN - www.suavinex.com - suavinex@suavinex.com - C.I.F. B03074093

	This appliance complies with the essential requirements established in the current EU directives.
	Caution
	Instructions for use
	In the EU this product should not be disposed of with other household waste in order to avoid harming the environment or human health. It must therefore be recycled correctly, promoting the sustainable reuse of material resources. To return the used appliance, please use the designated recycling systems or contact the establishment where you purchased the product. They will dispose of this product appropriately so that it can be safely recycled.
This product complies with EU Directive 2012/19/EU. The crossed-out wheeled bin symbol on the appliance indicates that the product must be disposed of separately from other household waste at the end of its useful life, either by taking it to a selective collection centre for electrical appliances or by returning it to the distributor when you purchase another similar appliance. The user is responsible for taking the appliance to a special waste disposal centre at the end of its useful life. If the unused appliance is correctly collected as separated waste, it can be recycled, processed and eliminated ecologically. This prevents it having a negative impact on the environment and health and helps ensure that the product materials are recycled. For more information on available waste disposal services, contact your waste disposal office or the shop where you bought the appliance. Google and Google Play are trademarks of Google LLC.	

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Google e Google Play são marcas comerciais da Google LLC.